




# PACKRACK PACKAGE TESTING FIXTURE

Testing of complete packages


# PACKRACK PACKAGE TESTING FIXTURE

PRODUCT BROCHURE

Package testing for  
bottles, caps, trays, cups,  
pouches, and more.

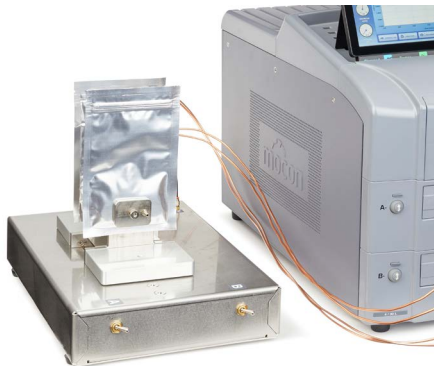

PackRack connected to an Instrument

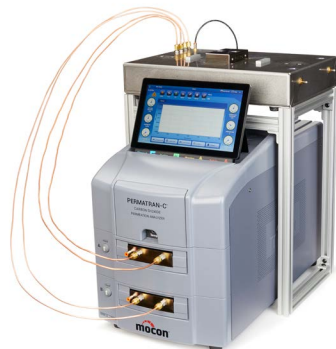

PackRack with space-saving stand option

## Complete Package Testing With Ease

The MOCON PackRack Package Testing Fixture expands the testing capability of MOCON film permeation analyzers from only 2 samples to 4 samples at a time for easy and accurate testing of complete packages or other 3-dimensional samples.

Revolutionary package testing cartridges simplify sample preparation and increase test result accuracy and repeatability. Some cartridge designs allow mounting of samples without the use of messy epoxy or glue. Others utilize inexpensive single use foils to mount the samples for repeatable results. These single use foils eliminate the messy and hazardous heating and scraping process used to clean the traditional mounting plates for re-use. The cartridges can be prepared in a separate designated workspace off-line to be ready for testing when needed. These same package testing cartridges are compatible with MOCON's other package testing permeation analyzers like the OX-TRAN 2/40, OX-TRAN 2/48 and the AQUATRAN 3/40.

## Features

- 2-cell fixture without any complicated electronics
- Enables film permeation analyzers to accurately test packages
- Less mess than traditional testing plates or custom test configurations
- Simplified package mounting without tools reduces labor and saves time
- Connects to most MOCON permeation analyzers with simple adaptor and tubing
- Can be used in environmental chambers

## Specifications

- Base model designed for benchtop.
- Available stand option allows mounting above an instrument to save space.
- Requires dedicated N<sub>2</sub> gas line for clamping function

We offer a variety of cartridge options designed to fit a wide range of package types.\*



Bottles, cups, caps, and more



Whole packages, syringes, tubing, and more



Large trays, cups, flexible pouches, bags and more

\*Note: analyzer and cartridges are sold separately.

**AMETEK**

**mocon**

MOCON, Inc.  
7500 Mendelssohn Ave. N  
Minneapolis, MN 55428 USA  
info.mocon@ametek.com  
www.ametekmocon.com